


# PIRANHHA

WALLED OFF NECROSIS TREATMENT

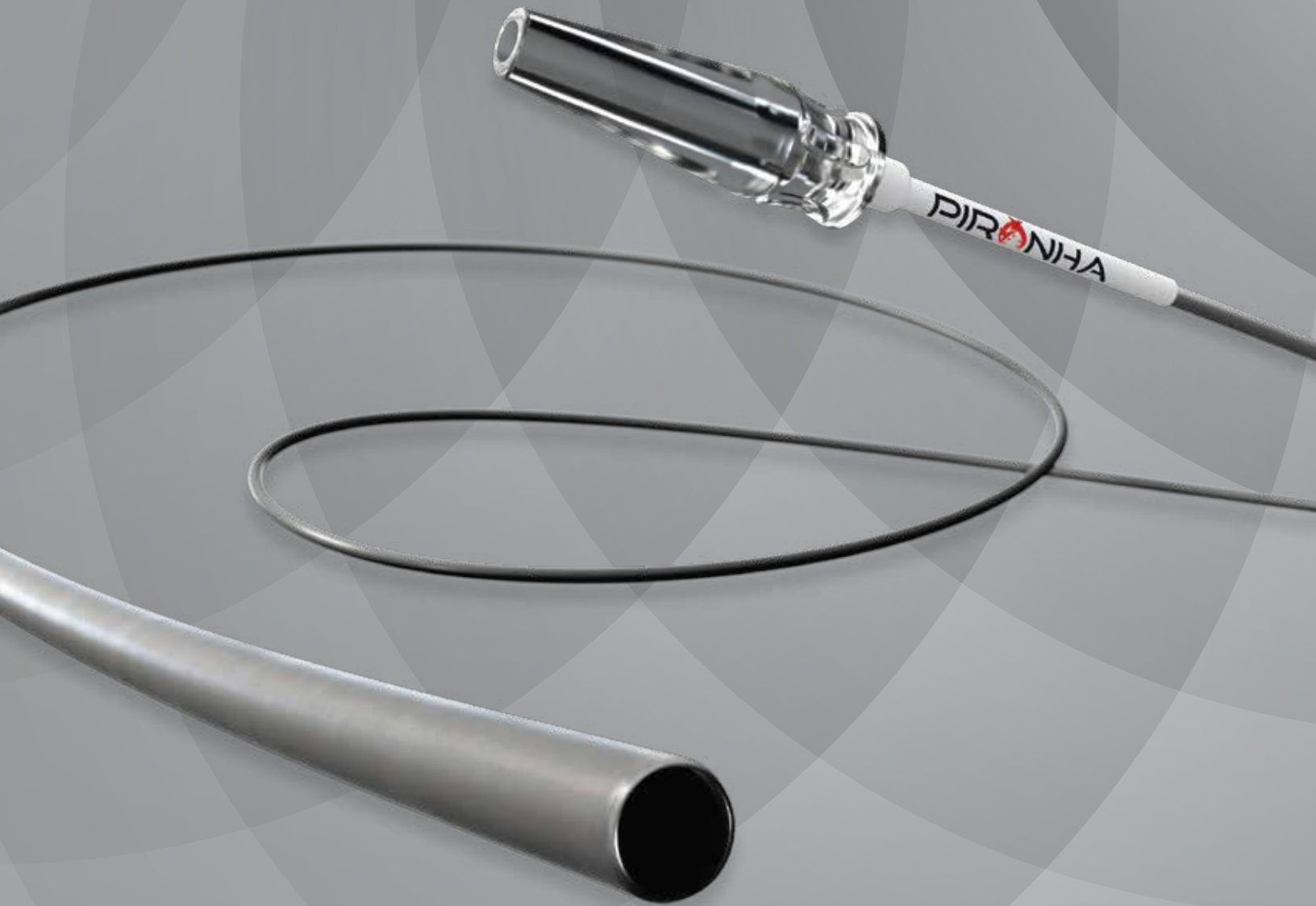

**ASSESS**  
pseudocyst

**REMOVE**  
fluid and tissue

**CLEAR**  
the field

[PiranhaMedical.com](http://PiranhaMedical.com)

The Piranha LC catheter is advanced through an endoscope



## DESIGN AND FUNCTIONALITY

Proximal 5" of shaft  
to add stiffness / pushability

Custom molded handle to  
interface with vacuum

Molded grips for allowing  
operator to hold at the  
handle or manipulate from  
the catheter shaft

Internal liner is Teflon coated  
providing a lubricious surface  
to maximize necrotic material  
suction

Custom tip with beveled edges and  
increased stiffness to aide in removal  
of necrotic material

